SJ-H series Digital Scale Operating Instructions

Cautions

- Avoid installing the scale in direct sunlight.
- Avoid over / shock load that could cause damage to the scale.
- Use the scale under conditions of -10°C ~ 40°C, less than 85% R.H.
- Do not mix the type of batteries. Replace the batteries all together.
- If the scale is not to be used for a long period of time, remove all batteries.
- Use only specified AC adapter for the scale (AC adapter is optional).

Specifications

SJ-1000H	1000 x 0.5 g / 2.2 x 0.001 lb / 35 x 0.02 oz / 2.2 lb x 0.1 oz
SJ-2000H	2000 x 1 g / 4.4 x 0.002 lb / 70 x 0.05 oz / 4.4 lb x 0.1 oz
SJ-5000H	5000 x 2 g / 11 x 0.005 lb / 176 x 0.1 oz / 11 lb x 0.1 oz
SJ-5001H	5000 x 1 g / 11 x 0.005 lb / 176 x 0.1 oz / 11 lb x 0.1 oz
SJ-12KH	12 kg x 5 g / 26 x 0.01 lb / 423 x 0.2 oz / 26 lb x 1 oz
Power	4 x "D" size batteries (not included) or AC adapter (option)
Battery life	Approximately 800 hours with alkaline cells at 20°C / 68°F
Pan size	230 (W) x 190 (D) mm / 9.1 (W) x 7.5 (D) in.
Dimensions	265 (W) x 250 (D) x 103 (H)mm / 10.4 (W) x 9.8 (D) x 4.1 (H) in.


Keys and Symbols

	key	Turns the scale power on/off.
	key	Sets the display to zero.
	key	Changes weighing units, (kg), lb, oz or lb-oz.
	Symbol	Indicates when the reading is stable.
	Symbol ZERO	Indicates when the scale zero is correct.
	Symbol NET	Indicates when net weight is displayed.
	Symbol	Battery remaining level, showed after power-on only.

Installing batteries / Connecting the AC adapter

Installing / Replacing Batteries

Slide the battery cover downward to open, and push the battery box to unhook and pull it out. Insert 4 batteries (D size) into the battery box, taking extreme care of the polarities.

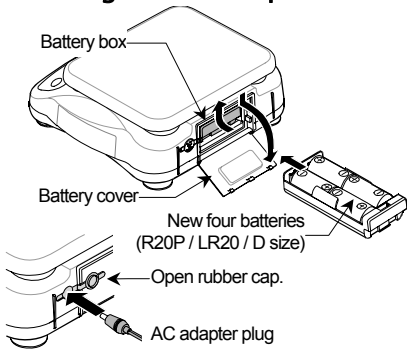AC adapter (optional)

Verify that the AC adapter is correct. Open the rubber cap and plug the AC adapter.

Weighing procedures

1. Turn on / off the scale

Press the key to turn on. The scale displays all segments for a few seconds and then shows zero. Press the key again to turn off.

2. Select the weighing unit

Press the key to select weighing unit. The scale will power up with the weighing unit that had been used when turned off last time.

3. Weighing

Verify the reading is zero. If not, press the key to zero the display. If you use a container for weighing, place an empty container on the weighing pan and press the key to zero the display. Place the objects to be weighed on the pan or in the container. When the reading becomes stable, the STABLE indicator is displayed.

Note

Automatic power-off function

If the scale is left ON and the STABLE is displayed, the automatic power-off function turns power off after approximately 5 minutes. To disable or enable this function, press the key while pressing the key. Then (disabled) or (enabled) will be displayed. Select the setting by pressing the key and press the key. The scale shows and returns to the weighing mode.

Error message

	Overload. Remove the objects from the pan.
	The scale is out of zero range or unstable when powered on.
	Low battery. Replace the batteries or use the AC adapter.

Calibration

Please visit the A&D website (<http://www.aandd.co.jp/>).

A&D Company, Limited

3-23-14 Higashi-Ikebukuro,
Toshima-ku, Tokyo 170-0013 JAPAN
<http://www.aandd.co.jp/>

A&D MERCURY PTY. LTD.

32 Dew Street, Thebarton, South
Australia 5031 AUSTRALIA
<http://www.andmercury.com.au/>

1WMPD4001461

A&D ENGINEERING, INC.

1756 Automation Parkway, San
Jose, California 95131 U.S.A.
<http://www.andonline.com/>

A&D KOREA Limited

8th Floor, Manhattan Bldg. 36-2
Yoido-dong, Youngdeungpo-ku,
Seoul, KOREA
<http://www.andk.co.kr/>

A&D INSTRUMENTS LTD.

Unit 24/26 Blacklands Way, Abingdon
Business Park, Abingdon, Oxon OX14 1DY
United Kingdom
<http://www.aandd-eu.net/>

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Notas

Calibration Instructions for the Type Approval Version of the SJ series

When is to Calibration Required?


Calibration may be required when the SJ series initially installed or if the scale is moved a substantial distance.

Calibration using a Calibration Weight

1. Entering calibration mode.

Remove the calibration switch cover at the bottom of the device and press the calibration switch while the device is in the Weighing Mode.

Then CAL will be displayed.

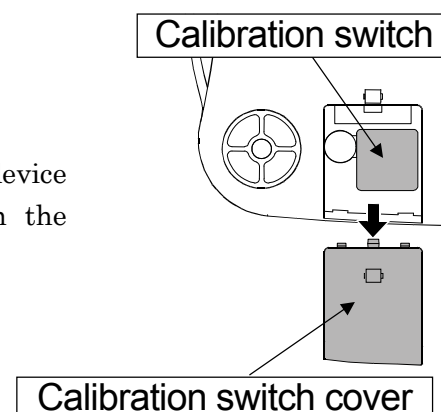2. Zero calibration

Press the [RE-ZERO] key to calibrate zero. CAL0 will be displayed.

Wait for the Stable Indicator to be displayed then press the [RE-ZERO] key.

CALF will be displayed after a few seconds. Press the [BLIND]* key to perform only a Zero Calibration. The scale will then automatically return to the Weighing Mode.

3. Span calibration

When CALF is displayed, place the weight at the center of the platform. Wait for the Stable Indicator to be displayed, then press the [RE-ZERO] key. The display will show End and the scale will automatically return to the Weighing Mode.

See "SPECIFICATIONS" about the calibration weight.

Calibration by Gravity Compensation

If the acceleration of gravity at your location is not 9.798 m/s² and you do not have calibration weights, the scale can be calibrated by compensating for the acceleration of gravity. (Refer to "The Value of Gravity at Various Locations".)

1. Setting a new acceleration value.

When CAL is displayed, press the [BLIND]* key and 9.798 will be displayed. Press the [RE-ZERO] key to increment the blinking digit and press the [BLIND]* key to move the blinking digit.

2. Storing the value in the memory.

While pressing the [BLIND]* key, press and hold the [RE-ZERO] key and release the [BLIND]* key. Then the display will show End and return to CAL. Turn the scale off to finish the procedure.

Remarks*: The [BLIND] key can be activated by pressing the area between the [ON/OFF] and [RE-ZERO] keys.

SPECIFICATIONS

MODEL	SJ-1000	SJ-2000	SJ-5000	SJ-12K
Capacity/Resolution	1000 g x 1 g	2000 g x 2 g	5000 g x 5 g	0 ~ 10 kg x 10 g 10 ~ 12 kg x 20 g
Maximum tare	1000 g	2000 g	5000 g	9.99 kg
Calibration weight	1000g ± 0.1g	2000g ± 0.2g	5000g ± 0.5g	10 kg ± 1g

The Value of Gravity at Various Locations

Amsterdam	9.813 m/s ²	Havana	9.788 m/s ²	Rio de Janeiro	9.788 m/s ²
Athens	9.807 m/s ²	Helsinki	9.819 m/s ²	Rome	9.803 m/s ²
Auckland NZ	9.799 m/s ²	Kuwait	9.793 m/s ²	San Francisco	9.800 m/s ²
Bangkok	9.783 m/s ²	Lisbon	9.801 m/s ²	Singapore	9.781 m/s ²
Birmingham	9.813 m/s ²	London (Greenwich)	9.812 m/s ²	Stockholm	9.818 m/s ²
Brussels	9.811 m/s ²	Los Angeles	9.796 m/s ²	Sydney	9.797 m/s ²
Buenos Aires	9.797 m/s ²	Madrid	9.800 m/s ²	Taichung	9.789 m/s ²
Calcutta	9.788 m/s ²	Manila	9.784 m/s ²	Vancouver, BC	9.809 m/s ²
Cape Town	9.796 m/s ²	Melbourne	9.800 m/s ²	Washington DC	9.801 m/s ²
Chicago	9.803 m/s ²	Mexico City	9.779 m/s ²	Wellington NZ	9.803 m/s ²
Copenhagen	9.815 m/s ²	Milan	9.806 m/s ²	Zurich	9.807 m/s ²
Cyprus	9.797 m/s ²	New York	9.802 m/s ²	Taiwan	9.788 m/s ²
Djakarta	9.781 m/s ²	Oslo	9.819 m/s ²	Taipei	9.790 m/s ²
Frankfurt	9.810 m/s ²	Ottawa	9.806 m/s ²	Tokyo	9.798 m/s ²
Glasgow	9.816 m/s ²	Paris	9.809 m/s ²		